


Vibrating Sample Magnetometer

High Field Magnet Laboratory

Radboud University

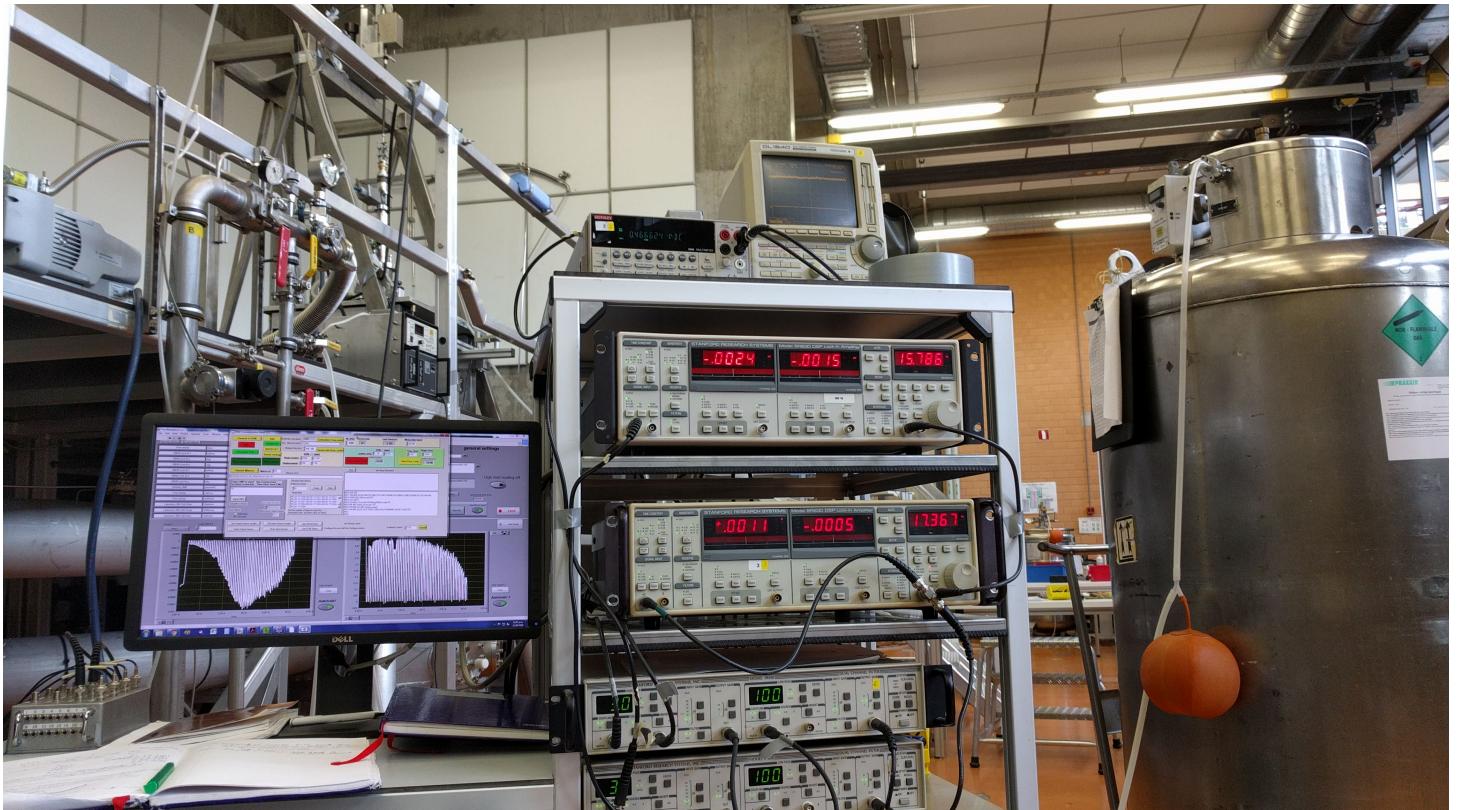

Andrew Ammerlaan

Supervisor: Prof. Dr. Uli Zeitler

November 7, 2018

Abstract

Vibrating sample magnetometry is a technique used to measure the magnetization of a sample. It utilizes the principle of electromagnetic induction; a magnetized sample will induce a voltage across a set of 4 pickup coils. This report contains both a in-depth analysis of the inner workings of the vibrating sample magnetometer (VSM) at HFML, and doubles as a manual for this VSM. Detailed schematics and illustrations with distance measurements are therefore included in the appendices. Finally, the VSM was used to measure the saturation magnetization of copper sulfate pentahydrate, using nickel as calibration.

Contents

1 Theoretical principle of the VSM	3
1.1 Magnetization and Faraday's law of induction	3
1.2 The response of a pickup coil to magnetized sample	3
1.3 Flattening the curve	4
1.4 Compensating for motion of the pickup coils	4
2 Setup	7
2.1 The drive motor	7
2.2 The pickup coils	7
2.3 The potentiometer	7
3 Calibration	8
3.1 Calibrating the potentiometer	8
3.2 Centering the sample	8
3.3 Calibrating with nickel	10
4 Saturation magnetization of CuSO₄·5 H₂O	11
5 Discussion	13
5.1 Different potentiometer power supplies	13
5.2 Magnetization at B is not the same as at -B?	13
5.3 Possibility for compensating for magnetization of empty sample holder	13
Acknowledgments	14
References	14
A Using the VSM motor controller program	15
B VSM measurement step-by-step	16
C Schematics and Pictures	18

1 Theoretical principle of the VSM

Materials in a magnetic field will react to this field, how they react exactly is an interesting field of study because not all materials will react in a trivial way. Some materials will exhibit phase transitions at certain magnetic field strengths, and some material's magnetization will have an interesting temperature dependence. By analyzing the magnetization of samples as a function of temperature and/or magnetic field strength, we can not only uncover interesting phenomena but also learn more about the sample material.

A vibrating sample magnetometer is capable of doing these magnetization measurements, using a set of four pickup coils. Two coils are positioned above the sample and two below. The sample is placed in a sample holder attached to a long stick. This sample stick is then screwed onto a connector piece, which in turn screws onto the top of the VSM. The length of the sample stick is adjusted such that the sample is positioned exactly between the pickup coils. The whole VSM is then inserted into a cryostat for temperature control, which in turn is inserted into a bitter magnet. This bitter magnet induces a magnetization on the sample. The connector piece on top of the VSM is attached to a motor, which sends a vibrating motion through the connector piece and causes the sample to vibrate in between the coils. The following sections explain the physical effect of the vibrating magnetized sample on a pickup coil, why 4 pickup coils are required, and why the magnetization can thus be measured accurately. Schematics and pictures of the VSM can be found in Appendix C.

1.1 Magnetization and Faraday's law of induction

Faraday's law describes the voltage that is induced in a closed loop due to a change in the magnetic field through this loop. It is given by:

$$V = -\frac{d\phi_B}{dt} \quad (1)$$

Where V is the induced voltage (also known as the 'electromotive force'), and ϕ_B is the magnetic flux through the loop which is given by:

$$\phi_B = \iint_S \vec{B}(\vec{r}, t) \cdot d\vec{A} \quad (2)$$

Where S represents the surface of the closed loop and B the magnetic field. By taking the magnetic field to be perpendicular to the surface of the closed loop, and homogeneous throughout the closed loop, Equation 2 simplifies to $\phi_B = BS$

Magnetization: The electrons in any material will react to an external magnetic field applied to this material, it will do this paramagnetically and/or diamagnetically.

Paramagnetic contributions are caused by the **magnetic moment of unpaired electrons** in the material that will align with the external magnetic field. This causes an internal magnetic field in the material which is parallel to the external field (positive magnetization).

Diamagnetism on the other hand is caused by the **Lorentz force on the electrons** in their orbitals causing them to create an internal magnetic field antiparallel to the external field (negative magnetization). All materials are diamagnetic to some extent, however because the electron orbitals are constrained by the Pauli principle these contributions are generally much smaller than paramagnetic contributions.

Once a material is magnetized it will have its own magnetic field, and will thus cause a magnetic flux through a closed loop. When the material moves with respect to the closed loop the flux will change, inducing a voltage in the closed loop. Because the induced voltage is directly related to the magnetization, it is possible to determine a sample's magnetization just by measuring the induced voltage.

1.2 The response of a pickup coil to magnetized sample

Because the inductive effect on just one closed loop is small, and because the alignment of one closed loop with respect to the magnetic field is prone to errors, a pickup coil is used instead. Since a coil is basically just a series of closed loops, we can multiply Equation 1 with N , the number of windings:

$$V(t) = -N \frac{d\phi_B}{dt} \quad (3)$$

This can be written as:

$$V(z, t) = -N \frac{d(BS)}{dz} \frac{dz}{dt} \quad (4)$$

Where S is a property of the sample proportional to its magnetization, and the z direction is parallel to the vibrating motion [4].

In an ideal situation the motion of the sample is perfectly sinusoidal, in other words $z(t) = A \cos(2\pi ft)$, and thus $dz/dt = -2\pi f A \sin(2\pi ft)$. If we further assume that there are no defects or other imperfections in the coil we can also use the following for its magnetic field:

$$B(z) = \frac{\mu_0 I}{2} (z - L) \left[\ln \left(r_2 + \sqrt{r_2^2 + (z - L)^2} \right) - \ln \left(r_1 + \sqrt{r_1^2 + (z - L)^2} \right) \right] + \frac{\mu_0 I}{2} z \left[\ln \left(r_1 + \sqrt{r_1^2 + z^2} \right) - \ln \left(r_2 + \sqrt{r_2^2 + z^2} \right) \right] \quad (5)$$

Where r_1 and r_2 are the inner and outer radius of the coil respectively, and L is the length of the coil [4].

Taking the derivative of Equation 5 and filling in dz/dt , we get:

$$V(z, t) = -NASI \mu_0 \pi f \sin(2\pi ft) \left[\ln \left(r_2 + \sqrt{r_2^2 + (z - L)^2} \right) - \ln \left(r_1 + \sqrt{r_1^2 + (z - L)^2} \right) + (z - L)^2 \left(\frac{1}{r_2 \sqrt{r_2^2 + (z - L)^2} + r_2^2 + (z - L)^2} - \frac{1}{r_1 \sqrt{r_1^2 + (z - L)^2} + r_1^2 + (z - L)^2} \right) + \ln \left(r_1 + \sqrt{r_1^2 + z^2} \right) - \ln \left(r_2 + \sqrt{r_2^2 + z^2} \right) + z^2 \left(\frac{1}{r_1 \sqrt{r_1^2 + z^2} + r_1^2 + z^2} - \frac{1}{r_2 \sqrt{r_2^2 + z^2} + r_2^2 + z^2} \right) \right] \quad (6)$$

From this it is immediately clear that the voltage signal from the pickup coil is also sinusoidal (90° phase shifted). In a experimental setup we will always measure the RMS value of this AC signal, therefore we take the RMS value of Equation 6:

$$V(z) = -NASI \mu_0 \frac{f \pi}{\sqrt{2}} \left[\ln \left(r_2 + \sqrt{r_2^2 + (z - L)^2} \right) - \ln \left(r_1 + \sqrt{r_1^2 + (z - L)^2} \right) + (z - L)^2 \left(\frac{1}{r_2 \sqrt{r_2^2 + (z - L)^2} + r_2^2 + (z - L)^2} - \frac{1}{r_1 \sqrt{r_1^2 + (z - L)^2} + r_1^2 + (z - L)^2} \right) + \ln \left(r_1 + \sqrt{r_1^2 + z^2} \right) - \ln \left(r_2 + \sqrt{r_2^2 + z^2} \right) + z^2 \left(\frac{1}{r_1 \sqrt{r_1^2 + z^2} + r_1^2 + z^2} - \frac{1}{r_2 \sqrt{r_2^2 + z^2} + r_2^2 + z^2} \right) \right] \quad (7)$$

We can check if this equation holds by making a plot of the response of 1 pickup coil as a function of the position at which the sample is vibrating (z), and fitting Equation 7 through this data, as shown in Figure 1.

1.3 Flattening the curve

The attentive reader will notice 2 things in Figure 1. First, the quantity on the y-axis: This is actually not $V(z)$ as in Equation 7 but $V(z)$ divided by a reference voltage that is directly proportional to the vibrating motion of the VSM's motor. The reason for this will be discussed in more detail in Section 2.3. Second, the plot in Figure 1 is not flat anywhere, which is a problem because this means that the sample will induce a different response throughout its vibrating motion. Because of this the signal is not going to be a perfect sine, since the signal will be higher closer to the coil and lower further away. For a accurate measurement we a nice sine wave as VSM signal, thus we need a flat area in the plot that is at least as big as twice the amplitude of the sample's motion.

This can be accomplished by measuring the combined responses of several pickup coils instead of just one. When a second coil is added on the other side of the sample and connected such that it responds oppositely (its response curve (Figure 1) is mirrored in the z-axis) the peaks of both coil start to merge, as demonstrated in Figure 2. By carefully choosing the separation between these 2 coils, we get a nice and flat (0.1% deviation [4]) area, as shown in Figure 3. (This is actually the combined response of not 2 but 4 pickup coils, see also Section 1.4)

1.4 Compensating for motion of the pickup coils

It is impossible to keep the VSM completely still inside the bitter magnet. There will always be unwanted vibrations caused by the various pumps, motors or water currents throughout the lab. These unwanted

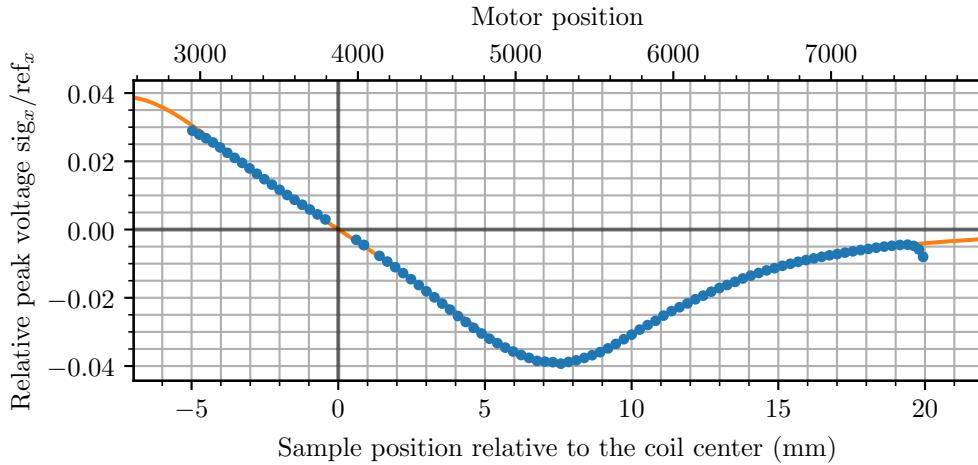


Figure 1: The response of one of the pickup coils (upper large, designated '1-2' in Appendix C on page 21) to the vibrational motion of a permanent magnet with a mass of 30.9 mg as a function of the position. Note that the coordinates on the lower axis have been shifted such that 0 responds to the coil center. The theoretical curve has been fitted through the data, Equation 7.

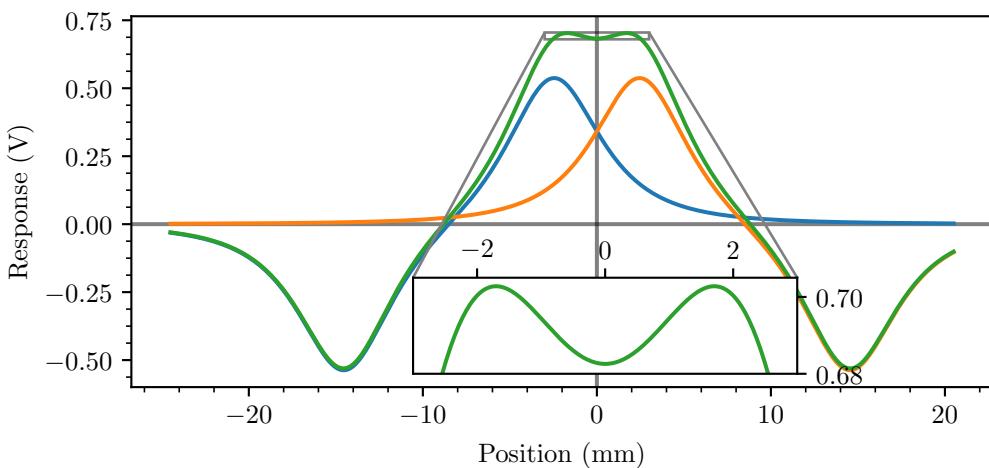


Figure 2: The superposition of the response of 2 coils, separated by a distance of 5 mm. Where the blue and orange curve represent the different coils, and the coils are connected such that the bottom of both coils are connected to each other. Because both coils are on opposite sides of the sample, this will cause both coils to respond with the same sign. In other words, one of the graphs is mirrored in the z-axis such that both peaks superimpose (green curve). Creating a approximately flat area on the peak of the superposition of both coils. For demonstration purposes the separation between the coils has been increased to 5 mm, the actual separation is 2.5 mm. Figure 3 contains the equivalent of the green curve with actual data points, but for a shorter z-axis range.

vibrations cause the pickup coils themselves to move inside the external magnetic field. While this field is quite homogeneous (1×10^{-3} cm DSV in cell 3 [2]) any such movement results in a small magnetic flux that will result in noise on the signal. A solution is to add an oppositely wound twin coil inside both pickup coils. This can compensate for any movements inside the external magnetic field, since these unwanted vibration will cause both this compensation coil and its parent coil to move in exactly the same way. For this compensation coil to cancel out these external vibrations the flux has to be exactly opposite to the flux that these vibrations cause in the parent coil. To accomplish this the compensating coil has to have a smaller radius, and thus needs to have a larger number of windings, such that the surface area of one winding multiplied by the number of windings is equal. When the flux through both

coils cancel each other out completely, there will be no overall effect of external vibrations on the signal when measuring across the combination of all coils.

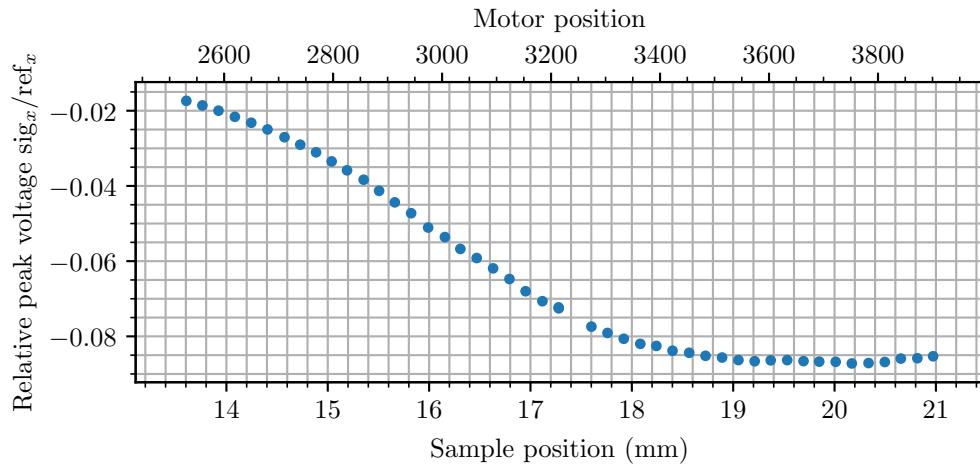


Figure 3: The response of the combination of all four pickup coils (designated '9-10' in Appendix C on page 21) to the vibrational motion of a nickel sample with a mass of 42.5 mg as a function of the position. The magnetization has been induced on the sample by an external magnetic field of 2 T in cell 3. The numbers on the lower axis are converted from the ones on the upper axis, this axis has its 0 when the motor driving the VSM is in the all the way up position.

Note that Figure 3 covers a shorter z-axis range than Figure 1, due to restrictions in the movement of the motor driving the VSM (more on this in the Section 2.1 and Section 2.3). If we were to continue the measurement beyond 21 mm, it will of course approximately follow the curve in Figure 2 and return to zero.

2 Setup

The whole VSM setup consists of 3 major components: the motor, the pickup coils and a potentiometer. This section will look into each component, and additional equipment needed to correctly analyze and record data from each component. A detailed schematic of the VSM setup can be found in Appendix C on page 21

2.1 The drive motor

The samples motion is driven by an electric motor, which is in turn controlled by a motor controller (LAC-25, Figure 13). This motor controller can receive commands from a computer using a VSM control program written by Hung van Luong. In this control program the user sets the desired amplitude and starting position, after which the computer will send the necessary commands to the motor. More detailed instructions on how to use this program can be found in Appendix A.

However if the motor were to be connected directly to the sample stick, the pumping up/down motion of the stick will cause undesired pressure changes inside the VSM. Therefore the motor does not drive the sample's motion directly instead it drives a connector piece on top of the VSM (Figure 12a). This connector piece seals the inside of the VSM and therefore prevents any pressure changes. The connector piece then passes on the vibrating motion of the motor to the sample stick (Figure 12b).

2.2 The pickup coils

When a sample is set to vibrate such that the sample is within the flat area (in Figure 3), we get a sinusoidal VSM signal from the combination of all 4 pickup coils, the frequency of which should be the same as the frequency of the sample's vibration. Previous research [3] has determined that a frequency of 18.1 Hz causes least noise in the VSM's signal, so this will be the frequency we will send to the motor. The VSM signal first goes through a Stanford Research 650 Dual-Channel filter, set to filter out everything below 4 Hz and above 80 Hz. Furthermore the filter also amplifies the signal, the amplification required depends on the amplitude of the original signal, and thus on the magnetization. The signal is then sent through a Stanford Research 830 Lock-In amplifier, which measures both the real and imaginary components of the signal and send these back to the computer. The signal can also be monitored visually on a oscilloscope, using the 'monitor out' of the lock-in. This is useful, because it enables us to check if the signal is truly sinusoidal and allows us to choose the correct amplification.

2.3 The potentiometer

Because of friction the actual motion of the motor is not necessarily equal to what was set by the user on the computer. This would not be a big problem if the friction were to be the same at each motor position. However there is a position dependent component to this friction, caused by the (springs of the) connector piece whose friction depends strongly on how much the spring is compressed (Section 2.1 and Figure 12a). As a result the amplitude when vibrating in a position that leaves the connector piece's springs more or less in their rest position will be significantly larger than when the connector piece's springs are pushed all the way down or up. This reduction in amplitude is directly visible in the VMS's signal (right part of Figure 5), which might lead one to conclude that the magnetization has decreased, when in fact it has not.

To record the actual motion of the motor (and thus the sample) we connect a potentiometer (Honeywell SLF01N1500B6A) to the motor. The potentiometer will return a certain fraction of a voltage applied to it depending on the position of the motor (and thus the sample). When the motor starts vibrating this fraction will change accordingly, thus creating a AC (vibration) signal on top of a DC (position) signal. To supply the potentiometer with power, I have made a box (Figure 4, and Appendix C on page 21) that uses 2 AA batteries to create a 3 V voltage across the potentiometer. The box will also separate both the AC and DC components from the signal returned by the potentiometer. The DC component is then sent to a Agilent 34410 multimeter, which sends it back to the computer. The AC component is sent through a Stanford Research 650 Dual-Channel filter and Stanford Research 830 Lock-In amplifier first, after which it can be recorded on the computer and displayed on the oscilloscope as well. We then use the potentiometer as a reference by dividing the VSM's signal by the potentiometer's signal. Because both signals will reflect the changed amplitude in the same way dividing them like this cancels out the reduction in amplitude, thus eliminating the effects of the position dependent friction, as has been done in Figure 1 and Figure 3.

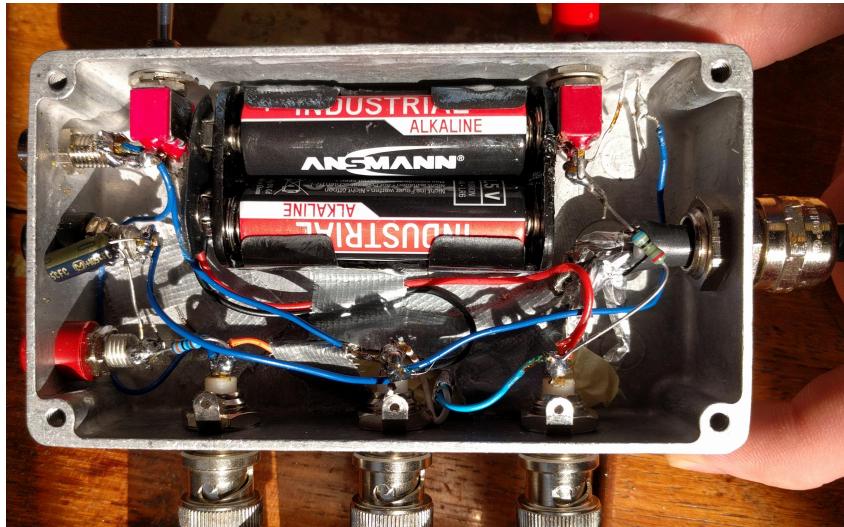


Figure 4: The box that supplies the potentiometer with power (2x1.5 V), and outputs the AC and DC component of the returned signal.

3 Calibration

Now the sample motion can be controlled, and the response signal from the pickup coils and the potentiometer can be measured. However, before any measurement can be done, the VSM first needs to be calibrated. First the potentiometer's DC component must be calibrated, in order to know the position of the sample at any given time. Next we need to find the position that leaves the sample exactly between the pickup coils, the flat area in a position loop plot. And finally the VSM's signal has to be calibrated with a sample of known magnetization.

3.1 Calibrating the potentiometer

To calibrate the potentiometer's DC component, simply set the VSM to vibrate at some position and read off the DC voltage on a multimeter or the computer. Do this for say 4 to 5 points, and make a linear fit through these points, this will yield parameters 'a' and 'b' ($position = a*voltage_{DC} + b$). These parameters will be used in the next section, to convert the recorded DC component back to a position and plot the VSM's signal as a function of this position (Figure 3). Note that this calibration needs to be redone every time the sample changes or the motor/connector piece/potentiometer is adjusted, because there is no way the (new) sample is going to be in exactly the same position after any adjustments.

3.2 Centering the sample

Next we must make sure that the tip of the VSM (coils + sample) is approximately in the center of the external magnetic field. This can be done by using the given specifications of the magnet [2], and the known lengths of the VSM parts (as found in Appendix C). The VSM can be moved up/down as needed with the screw shown in Figure 11. Furthermore to get a nice sinusoidal signal from the VSM it is even more important that the sample is set to vibrate exactly between the upper and lower pickup coils (the flat area in Figure 3). To accomplish this we need to measure the VSM signal as a function of the sample's position. For this purpose the VSM's control program has a 'position loop' option. This will set the VSM to vibrate at a set begin position for a set amount of time. After which it will move a set distance down, and vibrate for that same amount of time at that position, and so on.

When analyzing the VSM's signal on the lock-in we will see that most of the signal is imaginary since the coils impedance is mostly inductive. We can use the lock-in's "autophase" function to phase shift the signal such that almost all of the signal is real. We can then use the real part in our response versus position plots, since any change in the imaginary part is just noise.

Doing such a position loop for nickel will yield a graph as shown in Figure 5. In this figure we see that the resistive part (imaginary after "autophase") will fluctuate a bit as well. We can ignore this because we are only interested in the inductive part when looking for the flat area.

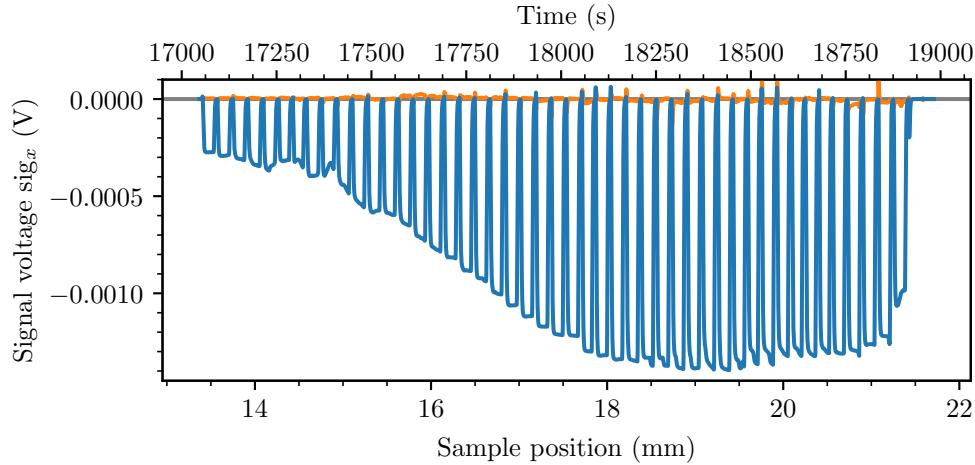


Figure 5: The response of the combination of four pickup coils (designated '9-10' in Appendix C on page 21) to the vibrational motion of a nickel sample with a mass of 42.5 mg as a function of the position. The magnetization has been induced on the sample by an external magnetic field of 2 T in cell 3. This figure contains the unprocessed data used to create Figure 3. The blue curve is the real part, while the orange curve represents the imaginary part.

If we then divide Figure 5 by the reference we get Figure 6. The advantage of dividing by the reference is immediately clear: Not only does the curve look smoother, it also compensates for the reduced amplitude of the VSM's motor (Section 2.3). This is particularly visible in the right of the plot, where in Figure 5 there is a significant reduction in amplitude. That is not visible at all in Figure 6.

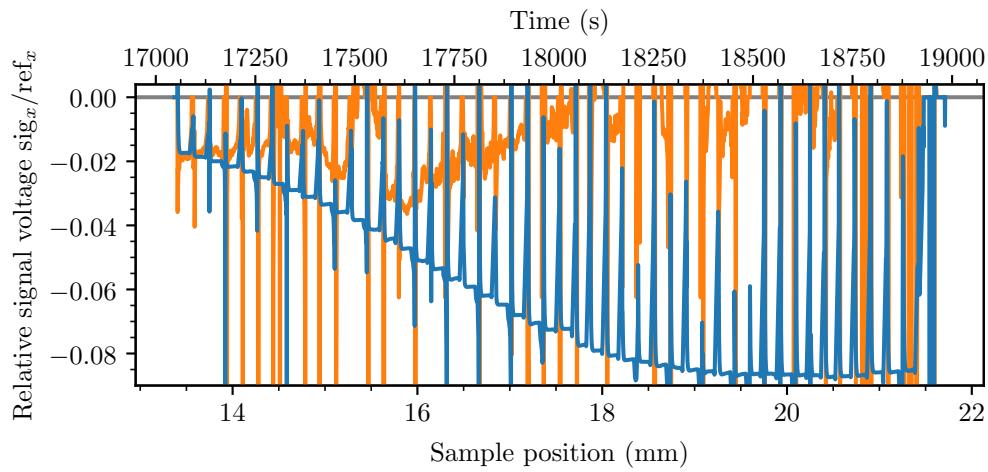


Figure 6: The response of the combination of four pickup coils (designated '9-10' in Appendix C on page 21) to the vibrational motion of a nickel sample with a mass of 42.5 mg as a function of the position. The magnetization has been induced on the sample by an external magnetic field of 2 T in cell 3. This figure is Figure 5 divided by the reference. The blue curve is the real part, while the orange curve represents the imaginary part.

However there is also a disadvantage to dividing by the reference, demonstrated by the spikes in Figure 6. These spikes are caused by points where the reference becomes almost 0 for a brief time, dividing by these near-zero numbers will cause very large positive or negative points in the plot. This is unavoidable since the motor has to stop before it can move to the next position. However, if we identify where each peak occurs in the data file before dividing by the reference, and use only those points after the division, we circumvent the divide-by-zero problem.

To automate this process I have made a simple python script, which can be found [here \[1\]](#). Which

takes data from any raw data file (e.g. Figure 5), locates each peak, divides by the reference, and then calculates the average and standard deviation of each peak. It then takes these data points and plots them against the position, converted from the DC component of the potentiometer's signal, and if able fits Equation 7 through the data (Figure 1 and Figure 3). More detailed instructions on how to use this can be found in the readme file, and on github [1]. Now that we have such a plot, we can easily identify the flat area and choose the sample's vibration position accordingly.

3.3 Calibrating with nickel

To measure the unknown magnetization of a sample, we must first calibrate the VSM's signal using a sample with a known magnetization, such as nickel. The saturation magnetization of nickel has been measured to be (58.57 ± 0.03) emu/g [5]. Therefore our sample (42.5 mg) should have a magnetization of (2.489 ± 0.002) emu

By setting the VSM to vibrate within the flat area and changing the strength of the external magnetic field, we obtain a plot as in Figure 7. Here we can see that nickel saturates at about 1 T, we can now calculate the average and standard deviation of all points from the saturation point onwards. We find that our saturated nickel sample $((2.489 \pm 0.002)$ emu) gives an absolute relative signal of 0.0848 ± 0.0002 . Now using the fact that 0 magnetization should give 0 signal, we can make a linear fit. Which gives a calibration constant of $c = (29.34 \pm 0.06)$ emu ($\text{magnetization} = c * \text{signal}$)

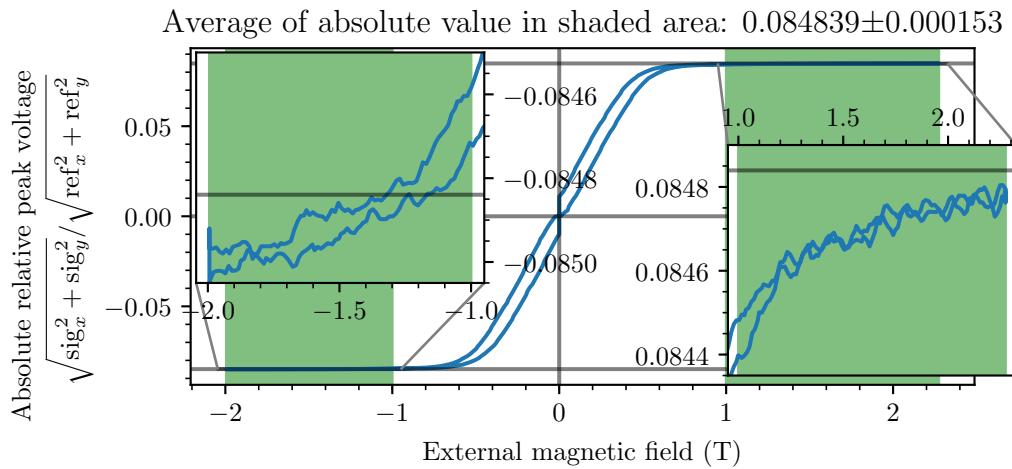


Figure 7: The response of the combination of four pickup coils (designated '9-10' in Appendix C on page 21) to the vibrational motion of a nickel sample with a mass of 42.5 mg. As a function of the applied external magnetic field in cell 3, at room temperature. The rate of change in the field was 20 mT/s, and the amplification on the signal and reference were both set to 30 dB.

4 Saturation magnetization of $\text{CuSO}_4 \cdot 5\text{H}_2\text{O}$

Now that we have calibrated the VSM, we can measure the unknown magnetization of any sample, such as copper sulfate pentahydrate. This sample has a very small magnetization, since it consists mostly of diamagnetic atoms, thus enabling us to test how sensitive the VSM is. In fact, at room temperature the magnetization of this material is so small, that it is near impossible to distinguish from an empty sample holder (Figure 8). Besides, at room temperature copper sulfate will not saturate within the range of the cell 3 bitter magnet (max 33 T).

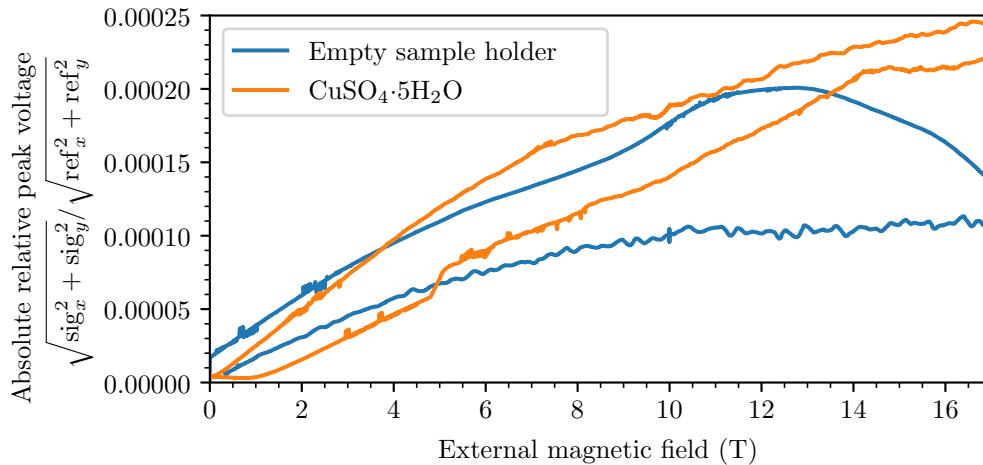


Figure 8: The response of the combination of four pickup coils (designated '9-10' in Appendix C on page 21) to the vibrational motion of a $\text{CuSO}_4 \cdot 5\text{H}_2\text{O}$ sample with a mass of 143 mg. As a function of the applied external magnetic field in cell 3, at room temperature. Compared to the same experiment with an empty sample holder.

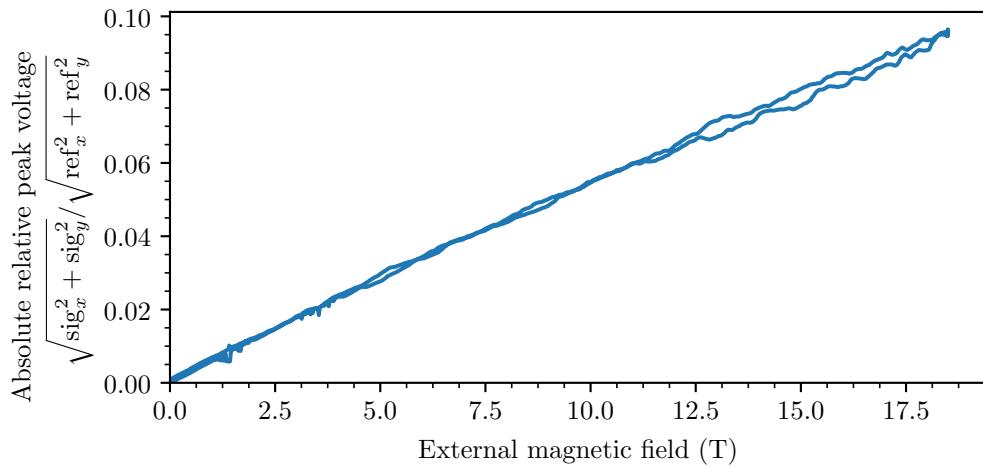


Figure 9: The response of the combination of four pickup coils (designated '9-10' in Appendix C on page 21) to the vibrational motion of a $\text{CuSO}_4 \cdot 5\text{H}_2\text{O}$ sample with a mass of 143 mg. As a function of the applied external magnetic field in cell 3, at 15 K. The rate of change in the field was 10 mT/s, and the amplification on the signal and reference were 40 dB and 20 dB respectively.

To measure the saturation magnetization anyway, we need to cool the sample down using cryostat 2 [2]. Because at low temperatures it is easier to induce magnetization, since there is less randomization due to thermal energy. The cryostat's thermometers (46682) and heater as well as the VSM's internal thermometer (x07377) are connected to a LakeShore 350 to control the cryostat. In Figure 9

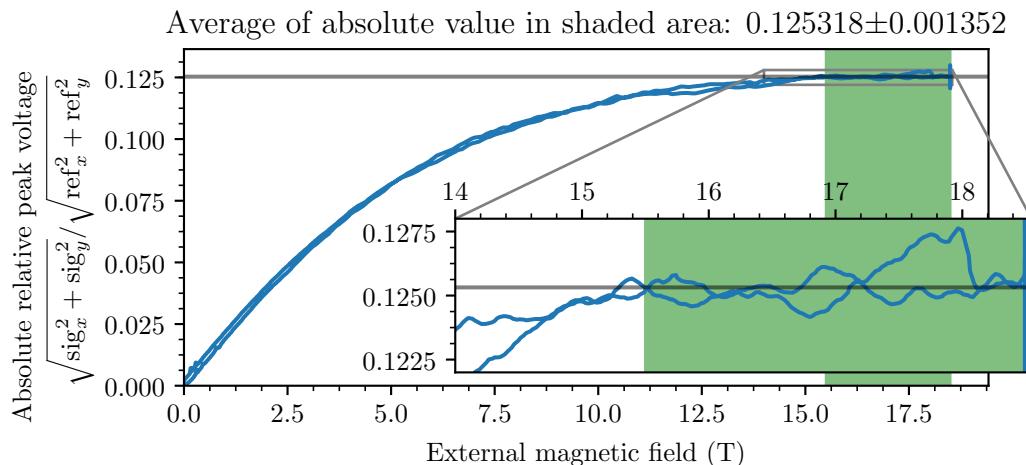


Figure 10: The response of the combination of four pickup coils (designated '9-10' in Appendix C on page 21) to the vibrational motion of a $\text{CuSO}_4 \cdot 5 \text{H}_2\text{O}$ sample with a mass of 143 mg. As a function of the applied external magnetic field in cell 3, at 4 K The rate of change in the field was 10 mT/s, and the amplification on the signal and reference were 40 dB and 20 dB respectively. The reference has been corrected to correspond to the reference used in Figure 7, using the slopes of the potentiometer calibration.

we see that copper sulfate does not yet saturate at 15 K, thus we need to go down further to 4 K (Figure 10) where saturation occurs at about 15.5 T. The calibration of the potentiometer was very different between the nickel and copper sulfate measurements ($y = 4287.347911 * x + 736.918625$ compared to $y = 8830.788615 * x + 238.4372817$, see section 5). This means that if we divide each measurement by their respective reference we cannot do a direct comparison between the two. However, since the calibration parameters for both measurements are known, we can easily compensate for this by using the ratio of the slopes (a_1/a_2) of both calibrations: After also compensating for amplification on the VSM's signal by the filter we can now divide the signal by the corrected reference and plot this (Figure 10). We see that at saturation we get a relative absolute signal of 0.125 ± 0.002 . Using the calibration parameter 'c' obtained in Section 3.3 we obtain a saturation magnetization of (3.67 ± 0.04) emu, which gives (25.7 ± 0.3) emu/g.

Sanity Check: In copper sulfate pentahydrate most of the magnetization is due to the Cu^{++} ions. These ions, being a spin $1/2$ system, have a spin magnetic moment of $2\sqrt{s(s+1)} = \sqrt{3} \mu_B$. Using the atomic mass given in table 1, we get a magnetization of 39 emu/g. Which is the same order of magnitude as the measured magnetization, they are however not equal. This is because we have only taken into account the magnetization of the copper ions, however the sulfate ions and water molecules will also have some diamagnetic contribution lowering the overall magnetization.

Element	Atomic mass (u)
Cu	63.546
S	32.065
9 O	9×15.9994
10 H	10×1.00794
$\text{CuSO}_4 \cdot 5 \text{H}_2\text{O}$	249.684

Table 1

5 Discussion

5.1 Different potentiometer power supplies

The large difference in potentiometer calibration in Section 4 arises from the fact that we have used different power supplies for the potentiometer throughout the duration of this report. Initially a Keithley 236 current source was used. However, this (very expensive) piece of equipment is a complete overkill for the purpose of supplying a voltage to a potentiometer. Therefore the plan was to make a box (Figure 4) that will supply this voltage, and which also includes a 1 Hz filter to separate the DC component.

We first tested the feasibility of this by using a simpler voltage source (Delta elektronika power supply E015-2), which was used during the copper sulfate experiments. However the voltage on this source was set higher than when using the Keithley 236 (during the nickel measurement). This lead to a significant difference between the 'a' parameters of both measurements.

The problem has been worked around in Section 4. However, using the calibration parameters in the calculation means that we need to take a look at the errors of these parameters. The parameters were obtained using linear regression, this also yielded the following coefficients of determination for nickel and copper sulfate respectively:

- 0.9999985027
- 0.9998095329

Thus the error in the linear fit is relatively small.

Note that when using the box as power supply for the potentiometer, it is still required to correct the reference if there is a large time delay between VSM calibration and measurement. This is due to the battery slowly depleting overtime, however if one records the current battery voltage for each measurement one can easily determine this correction factor.

5.2 Magnetization at B is not the same as at -B?

The attentive reader might have noticed that in Figure 7 at saturation the signal for positive B is always below the average-line, while for negative B it is above. This appears to suggest that the saturation magnetization of the nickel sample is larger for negative magnetic field and smaller for positive magnetic field. This is of course impossible, something else must be going on!

A possible explanation for this aberration would be a constant negative DC offset in the VSM's signal. However the signal goes through a 4 Hz high pass so any DC offsets should have been filtered out.

Another potential explanation is that there is some small contribution to the overall magnetization by the sample holder or impurities left behind by previous samples in the holder (we know the holder's magnetization is non zero, Figure 8). However these components would have to have a magnetization that is somehow not the same for different directions of the magnetic field.

A satisfying explanation to this strange behavior has yet to be found. However by including both the contributions at negative and positive field in the average and standard deviation we include the difference between the two in the error.

5.3 Possibility for compensating for magnetization of empty sample holder

In determining the calibration constant (Section 3.3) we assumed that 0 sample magnetization corresponds to 0 signal. However if we look at Figure 8 we see that this is not entirely true, as the sample holder (and impurities left behind by previous samples) also have some small (diamagnetic) contributions to the overall magnetization.

A possibility to compensate for this exists in using 2 different samples (e.g. Nickel and Iron) with known magnetization for calibration and linear fitting those ($y = ax + b$ instead of $y = ax$). This removes the assumption that the calibration-line should go through the origin ($b=0$), and thus allows for some signal to exist at 0 field without it interfering with the measurement of some unknown magnetization later on. Note that the empty sample holder leads to a signal of the order of 0.0002 (Figure 8) which is within the error (standard deviation: 0.002) of the copper sulfate measurement (Figure 10).

Acknowledgments

None of this would have been possible without the help of Prof. Dr. Uli Zeitler, who gave me the opportunity to work with the VSM at HFML. And was always there to answer my many questions, and to explain how all the things worked; the VSM, the lock-in amplifiers, the cryostat, the helium pump, and of course the bitter magnets. My sincerest thanks.

My thanks also to Hung van Luong, who taught me how to use the VSM control program and also showed my how all the devices (lock-ins, filters) should be connected to each other. And helped me whenever I had problems with the motor or the computer.

Furthermore, my thanks go out to Lijnis Nelemans, Michel Peters and Edwin van Leverink, who helped me create the potentiometer's box and provided me with technical advice and equipment.

I also would like to thank the rest of the amazing staff and students at HFML, who have been very kind and helpful in showing me around and helping me find all the things in the lab.

Thank you

References

- [1] A.A.L.N. Ammerlaan. *HFML-VSM pickupcoil profile plotter*. URL: https://github.com/AndrewAmmerlaan/HFML-VSM_CoilProfilePlotter.
- [2] *HFML magnet specifications*. URL: <https://www.ru.nl/hfml/facility/experimental/magnets/>.
- [3] O. Young. logbook.
- [4] E. Kampert. “Magnetic properties of organometallic compounds in high magnetic fields”. PhD thesis. Radboud University, 2012.
- [5] H. Danan, A. Herr, and A. J. P. Meyer. “New Determinations of the Saturation Magnetization of Nickel and Iron”. In: *Journal of Applied Physics* 39.2 (2008), pp. 669–670.

A Using the VSM motor controller program

To start motor:

- **Warning:** If you disconnect the motor's powersupply, the motor will have to be reinitialized ("Init"). This means that the motor will have to be disconnected from the VSM, which usually means that both the potentiometer's calibration and the VSM's center will change. Hence a full recalibration is required.
- Press "Connect to VSM", usually the VSM resides on COMM 3, and the default settings are fine
- Press "Ok"
- If not already on, press "Motor on"
- If this is first time the motor is used after it has been plugged in: press "Init". The motor will move as far up as it can go, and use that as the zero position. For this to work, the motor has to be disconnected from the rest of the VSM, and the potentiometer must be moved all the way up.
- Set the desired frequency, displacement (2x amplitude), and begin position (the motor will vibrate between 'begin position' and 'begin position' + 'displacement').
- Click "Calculate freq code", this will convert the frequency and displacement in a single code that the motor can understand.
- Click "Send calc freq code", which will send the frequency code to the motor.
- Click "Send string", this sends the string that contains the 'begin position' parameter to the motor.
- Click "Start motor", to stop the motor click "Stop motor".

Starting a position loop:

- Complete the above, but do not yet start the motor.
- Set the desired step size, number of steps and step duration.
- Click "Start position loop"
- When the loop has been completed, the motor will continue to vibrate at the last position, to stop the motor click "Motor on" immediately followed by "Motor off"

Troubleshooting:

- Motor does not respond:
The motor stores commands in finite volatile memory, when it is full the motor cannot receive any further commands from the controller program. To clear the memory, press the buttons "Flush input queue", "Flush output queue", and "Reset macro".
- Motor does not start:
If the motor does not start after it has been disconnected from the power, press "Init" to reinitialize the motor. Because the motor uses volatile memory, it will forget everything once the power has been disconnected.
- Motor still does not start:
Try to start the motor without anything attached to it. If it works now, it is because there was too much friction. Try adjusting the potentiometer, or turn the connector piece a bit.
- The motor stops halfway a position loop:
Beyond this point the motor encounters too much friction to continue, sometimes the friction can be reduced by turning the connector piece a bit for optimal alignment.

B VSM measurement step-by-step

Initial setup

- Put cryostat in magnet, and put VSM inside cryostat.
- Adjust VSM height (Figure 11) such that the pickup coils are in the center of the external magnet (Measure the distance from top of cryostat to the center of the magnet, and use schematics in appendix C).
- Connect lock-ins/amplifiers etc. as shown schematically on page 21.
- Connect cryostat to the lakeshore 350 and to the helium pump. Connect the VSM's thermometer to the lakeshore as well.
- Set the calibration curves the lakeshore should use for the thermometers (VSM:x07377, cryostat2:46682).
- Vacuum pump the cryostat's vacuum chamber.
- Setup measurement program to record signal, AC-reference, DC-reference, temperature and magnetic field
- Move potentiometer all the way up.
- Turn the motor on, and initialize (Appendix A), if not already initialized.
- Turn the motor off.

Putting a sample in the VSM

- Put calibration sample in sample holder, and screw it onto the sample stick (Figure 12b).
- Screw the sample stick onto the connector piece (Figure 12a), and put the stick in the VSM.
- Connect the motor to the connector piece.
- Adjust the connector piece such that the motor can move freely (Try screwing the connector piece tighter/looser), when done correctly the motor will fall down when pushed up and released.
- Move potentiometer such that the motor can make use of the full range of the connector piece.
- Turn motor on
- Find the smallest and largest possible motor position, do this by choosing some position and visually checking the connector piece.
- Calibrate the potentiometer, by setting the motor to a couple of different begin positions and writing down the corresponding DC component of the reference, and making a linear fit.

Cooling sample

- Fill cryostat with liquid nitrogen followed by liquid helium, depending on desired temperature.
- If required set desired temperature in the lakeshore and enabled the heater.
- Once desired temperature has been reached turn on the motor and the bitter magnet.

Doing a measurement

- Set the motor to vibrate at some position, and set the bitter magnet to some field where there is a clear signal, adjust signal/reference amplification if required. Note that for calibration samples that already have a magnetization the bitter magnet can stay off for the time being.
- Use the lock-in's "autophase" function
- Do a position loop from the smallest to the largest position. To save time it might be a good idea to first do a rough position loop with just a few steps, and finding where the center is approximately. Followed by a second more detailed position loop, with more steps, using an educated guess as begin position. Use coilprofileplotter.py to help find the center, detailed instructions in the readme file [1].
- Set the motor's begin position and displacement such that the sample will vibrate exactly in between the pickup coils.
- write down battery voltage.
- Measure as a function of the magnetic field.

Replacing sample

- Turn off the motor, and disconnect it from the connector piece.
- Create a over-pressure in the cryostat, unscrew the connector piece and remove the sample stick, then put the connector piece back on the VSM.
- Replace sample in sample holder.
- Unscrew the connector piece and reattach the sample stick, screw the connector piece back on the VSM and reconnect the motor.
- Wait for sample to reach desired temperature.
- Re-calibrate the potentiometer, and do a new position loop.

Calibrating VSM

- Steps "Putting a sample in the VSM" and "Cooling sample" for calibration sample.
- Steps "Doing a measurement".
- For better calibration: Steps "Replacing Sample" for second calibration sample.
- For better calibration: Steps "Doing a measurement".
- If necessary correct second reference using the battery voltages for the different measurements.
- Identify the saturation point(s) and find the relative absolute signal(s) at saturation.
- Use the known magnetization of the calibration sample(s) to obtain a calibration ' $c*x$ ' or ' $c*x+b$ '.

Measuring sample with unknown magnetization

- Steps "Replacing sample" for target sample.
- Steps "Doing a measurement"
- Identify saturation point and find the relative absolute signal at saturation, make sure to correct the reference first using the ratio of the battery voltages (if changed).
- Use calibration to find magnetization and magnetization per gram.

C Schematics and Pictures

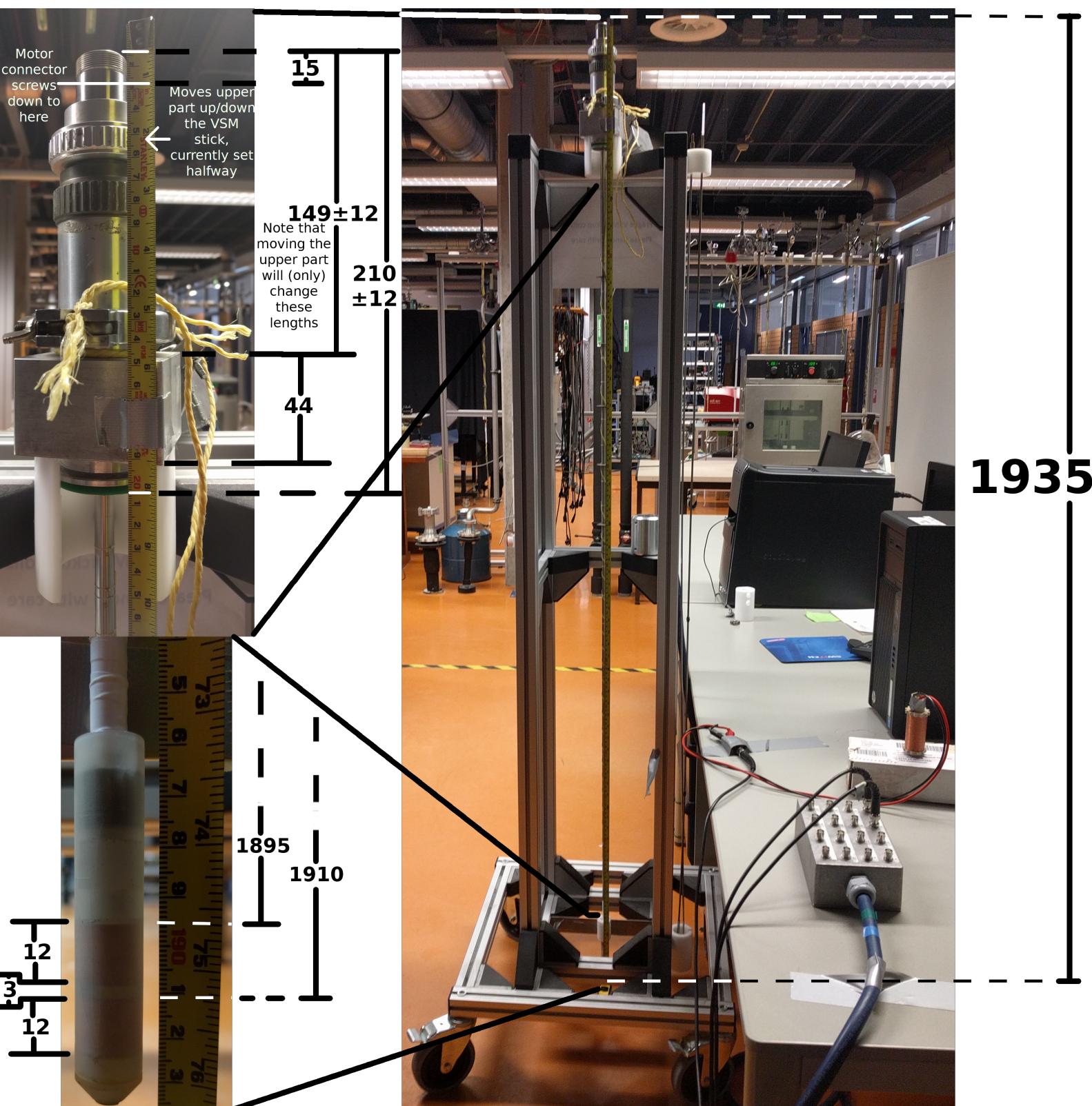


Figure 11

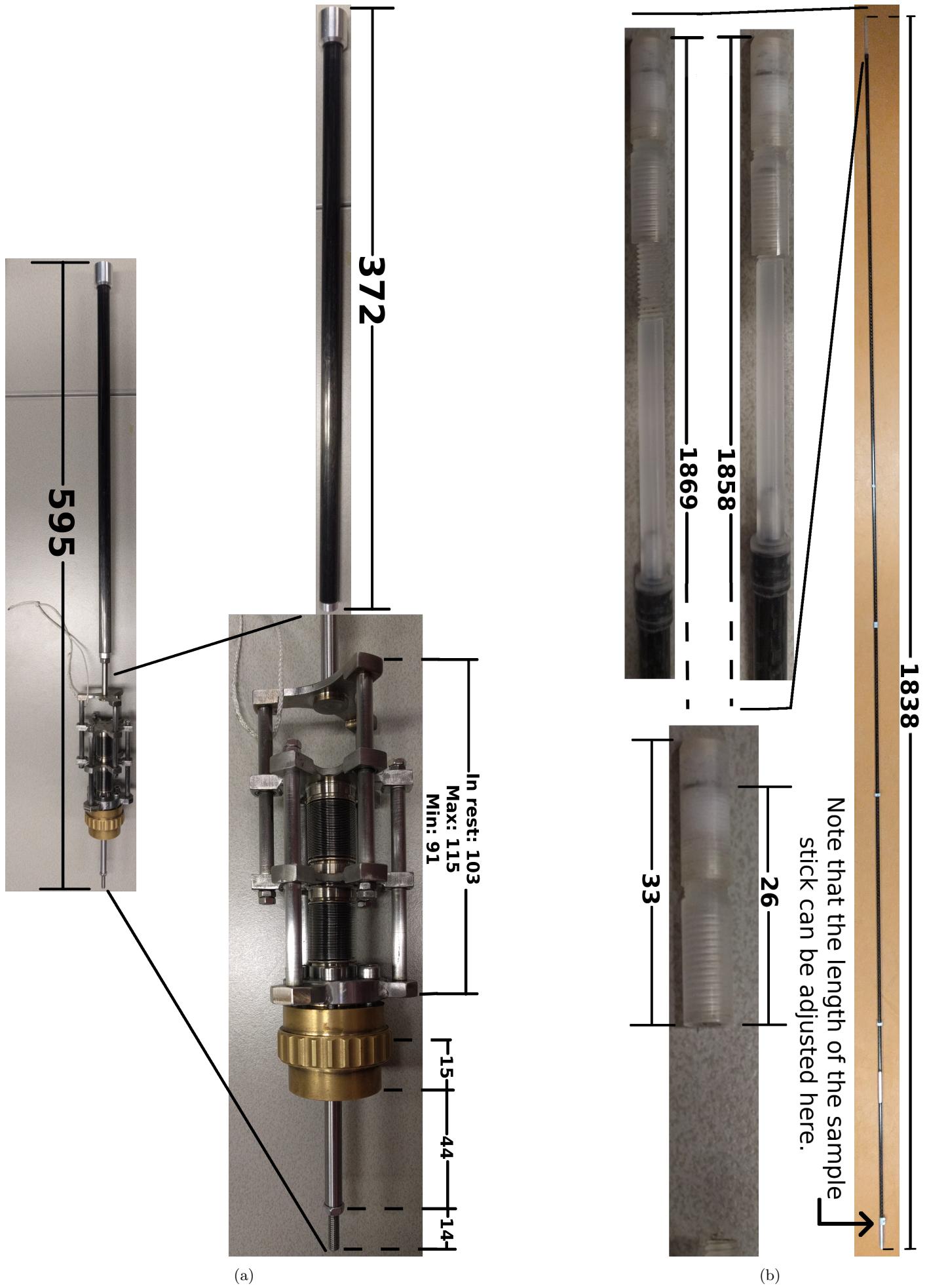


Figure 12

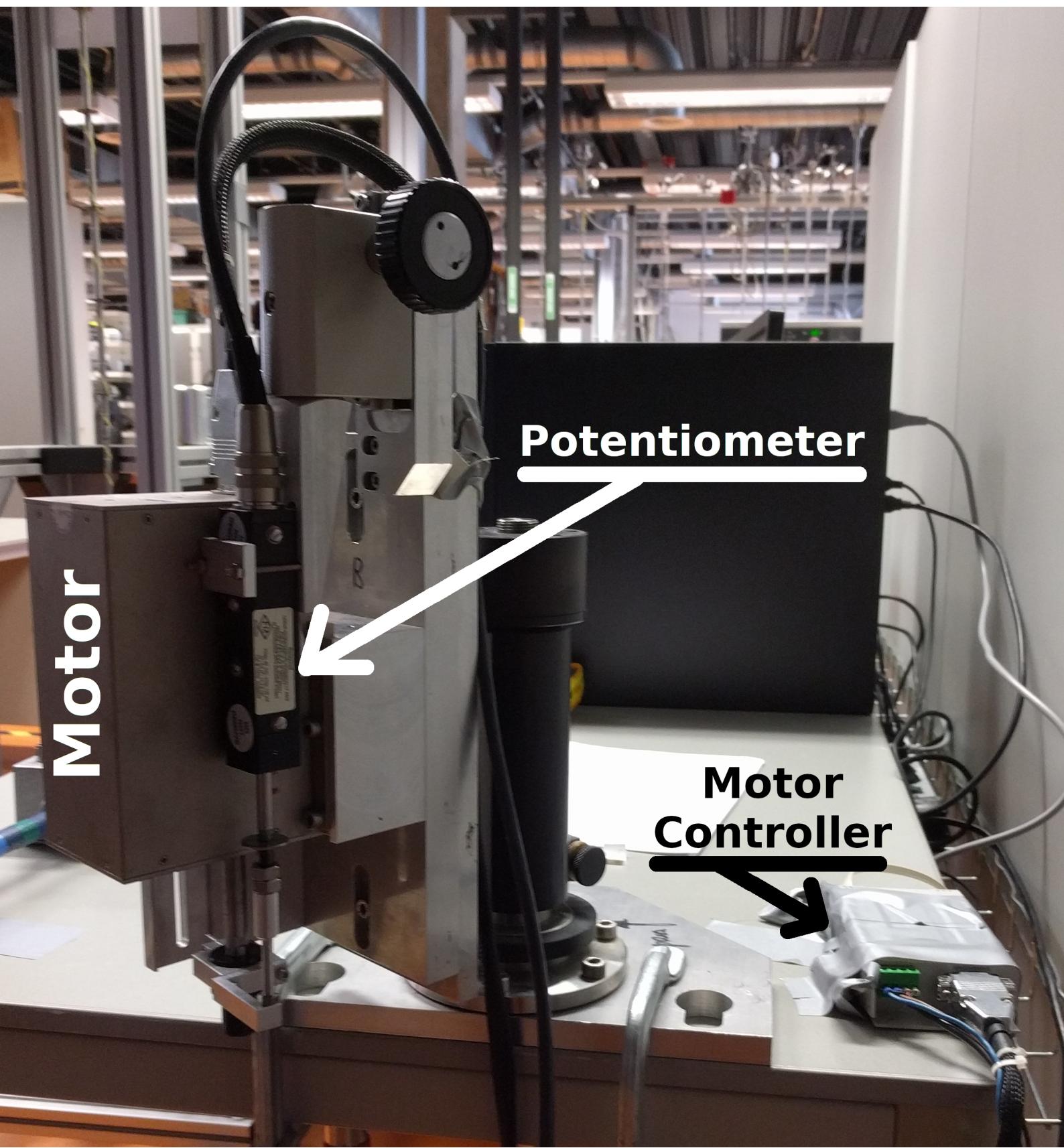


Figure 13

